Datasheet

Version: 1.0.0 Revision date: 23 Feb 2025



Pig Proprotein Convertase Subtilisin/Kexin Type 9 (PCSK9) ELISA Kit

Catalogue No.:abx353475

Pig Proprotein Convertase Subtilisin/Kexin Type 9 ELISA Kit is an ELISA Kit against Proprotein Convertase Subtilisin/Kexin Type 9.

Target: Proprotein Convertase Subtilisin/Kexin Type 9 (PCSK9)

Reactivity: Pig

Tested Applications: ELISA

Recommended dilutions: Optimal dilutions/concentrations should be determined by the end user.

Storage: Shipped at 4 °C. Upon receipt, store the kit according to the storage instruction in the kit's

manual.

Validity: The validity for this kit is at least 6 months. Up to 12 months validity can be provided on request.

Stability: The stability of the kit is determined by the rate of activity loss. The loss rate is less than 5%

within the expiration date under appropriate storage conditions. To minimize performance fluctuations, operation procedures and lab conditions should be strictly controlled. It is also

strongly suggested that the whole assay is performed by the same user throughout.

Test Range: 0.78 ng/ml - 50 ng/ml

Sensitivity: 0.47 ng/ml

Standard Form: Lyophilized

Detection Method: Colorimetric

Assay Type: Sandwich

Assay Data: Quantitative

Sample Type: Serum, plasma and other biological fluids.

Datasheet

Version: 1.0.0 Revision date: 23 Feb 2025



Kit Components:

The kit components listed are for reference only. The product manual may differ slightly. The product should be used as stated on the product manual included and delivered together with the product.

- Pre-coated 96-Well Microplate
- Standard
- · Standard Diluent Buffer
- Wash Buffer
- · Detection Reagent A
- · Detection Reagent B
- Diluent A
- Diluent B
- TMB Substrate
- Stop Solution
- Plate Sealer

Material Required But Not • 37°C incubator


Provided:

- Multi and single channel pipettes and sterile pipette tips
- · Squirt bottle or automated microplate washer
- 1.5 ml tubes
- Distilled water
- · Absorbent filter papers
- 100 ml and 1 liter graduated cylinders
- Microplate reader (wavelength: 450 nm)
- ELISA Shaker

Assay Procedure:

This procedure is provided for reference only. The product manual may differ slightly. The product should be used as stated on the product manual included and delivered together with the product.

- 1) Set standard, test samples and control wells.
- 2) Aliquot 100 µl of diluted standard into the standard wells.
- 3) Aliquot 100 µl of Standard Diluent buffer into control (zero) well.
- 5) Aliquot 100 µl of Detection Reagent A to each well. Incubate for 1 hr at 37 °C.
- 6) Wash 3 times.
- 7) Aliquot 100 µl of Detection Reagent B to each well. Incubate for 90 mins at 37 °C.
- 9) Aliquot 90 µl of TMB Substrate to each well. Incubate for 10-20 mins at 37 °C.
- 10) Aliquot 50 µl of Stop Solution.
- 11) Measure the OD at 450 nm.

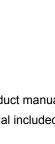



Datasheet

Version: 1.0.0 Revision date: 23 Feb 2025



Assay Precision:

Intra-assay Precision (Precision within an assay): 3 samples with low, medium and high levels of Proprotein Convertase Subtilisin/Kexin Type 9 were tested 20 times on one plate, respectively. Inter-assay Precision (Precision between assays): 3 samples with low, medium and high levels of Proprotein Convertase Subtilisin/Kexin Type 9 were tested on 3 different plates, 8 replicates in each plate.

CV (%) = (Standard Deviation / mean) × 100

Intra-Assay: CV<10% Inter-Assay: CV<10%

Note: This product is for research use only.

The range and sensitivity is subject to change. Please contact us for the latest product information. For accurate results, sample concentrations must be diluted to mid-range of the kit. If you require a specific range, please contact us in advance or write your request in your order comments.

Please note that our ELISA and CLIA kits are optimised for detection of native samples, rather than recombinant proteins. We are unable to guarantee detection of recombinant proteins, as they may have different sequences or tertiary structures to the native protein.



Website: www.abbexa.com · Email: info@abbexa.com